

# CDN Resource Laboratories Ltd.

#2, 20148 – 102<sup>nd</sup> Avenue, Langley, B.C., Canada, V1M 4B4, 604-882-8422, Fax: 604-882-8466 (www.cdnlabs.com)

## REFERENCE MATERIAL: CDN-GR-2

Recommended value and the "Between Laboratory" two standard deviations

<i>Graphitic carbon concentration</i>	<i>1.93 ± 0.21 % C</i>	<i>Certified value</i>
---------------------------------------	------------------------	------------------------

**PREPARED BY:** CDN Resource Laboratories Ltd.  
**CERTIFIED BY:** Duncan Sanderson, B.Sc., Licensed Assayer of British Columbia  
**INDEPENDENT GEOCHEMIST:** Dr. Barry Smee., Ph.D., P. Geo.  
**DATE OF CERTIFICATION:** March 20, 2014

### ORIGIN OF REFERENCE MATERIAL:

Standard CDN-GR-2 was prepared using ore supplied by Noram Ventures Inc. from their Kokanee Graphite property near Crawford Bay on Kootenay Lake in southeastern British Columbia. Metamorphic graphite is hosted in pelitic and epiclastic sedimentary rocks of the Paleozoic Index Formation of the Lardeau Group within the Kootenay Arc Metamorphic Complex. These rocks are amphibolite facies quartzites and quartz-muscovite-biotite ± garnet and sillimanite schists.

### METHOD OF PREPARATION:

Reject ore material was dried, crushed, pulverized and then passed through a 270 mesh screen. The +270 material was discarded. The -270 material was mixed for 5 days in a double-cone blender. Splits were taken and sent to 8 commercial laboratories for round robin assaying. Round robin results are displayed below:

### APPROXIMATE CHEMICAL COMPOSITION (by whole rock analysis):

	Percent		Percent
SiO <sub>2</sub>	57.5	Na <sub>2</sub> O	0.2
Al <sub>2</sub> O <sub>3</sub>	7.9	MgO	1.9
Fe <sub>2</sub> O <sub>3</sub>	3.3	K <sub>2</sub> O	2.9
CaO	11.6	TiO <sub>2</sub>	0.3
MnO	<0.1	LOI	10.9
Total S	1.0		

### ASSAY RESULTS

	Lab 1	Lab 2	Lab 3	Lab 4	Lab 5	Lab 6	Lab 7	Lab 8
Sample	Grph C %	Grph C %	Grph C %	Grph C %	Grph C %	Grph C %	Grph C %	Grph C %
GR-2-1	1.94	2.03	1.99	2.02	1.65	1.70	1.91	1.95
GR-2-2	1.91	2.01	1.98	2.03	1.70	1.90	1.89	2.06
GR-2-3	1.91	2.01	1.99	2.06	1.80	2.20	1.92	2.00
GR-2-4	1.96	2.04	1.98	2.05	1.85	1.90	1.96	1.99
GR-2-5	1.97	2.04	1.95	2.05	1.80	1.90	1.94	2.00
GR-2-6	2.00	2.06	1.96	2.05	1.85	1.80	1.92	2.03
GR-2-7	1.90	2.04	2.00	2.04	1.85	1.80	1.93	2.04
GR-2-8	1.79	2.05	1.98	2.05	1.90	1.90	1.86	1.99
GR-2-9	1.94	2.01	1.98	2.03	1.80	2.00	1.92	2.00
GR-2-10	1.98	2.05	1.93	2.02	1.85	2.00	1.93	1.97
Mean	1.93	2.03	1.97	2.04	1.81	1.91	1.92	2.00
Std. Dev'n	0.0591	0.0184	0.0212	0.0141	0.0762	0.1370	0.0274	0.0327
%RSD	3.06	0.90	1.07	0.69	4.22	7.17	1.43	1.63

**REFERENCE MATERIAL: CDN-GR-2**

**Statistical Procedures:**

The final limits were calculated after first determining if all data was compatible within a spread normally expected for similar analytical methods done by reputable laboratories. Data from any one laboratory was removed from further calculations when the mean of all analyses from that laboratory failed a t test of the global means of the other laboratories. The mean and standard deviation were calculated using all remaining data. Any analysis that fell outside of the mean  $\pm 2$  standard deviations was removed from the ensuing data base. The mean and standard deviations were again calculated using the remaining data. This method is different from that used by Government agencies in that the actual "between-laboratory" standard deviation is used in the calculations. This produces upper and lower limits that reflect actual individual analyses rather than a grouped set of analyses. The limits can therefore be used to monitor accuracy from individual analyses, unlike the Confidence Limits published on other standards.

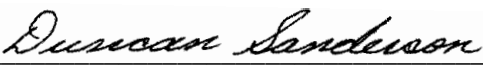
Participating Laboratories: (not in same order as table of assays)

Acme Analytical Laboratories Ltd., Vancouver, BC, Canada  
Activation Laboratories, Ancaster, Ontario, Canada  
ALS Canada, North Vancouver, B.C., Canada  
ALS Brisbane, Australia  
Amdel, Australia  
Intertek – Genalysis, Perth, Australia  
SGS, Vancouver, BC, Canada  
SGS – Lakefield, Ontario, Canada


Legal Notice:

This certificate and the reference material described in it have been prepared with due care and attention. However CDN Resource Laboratories Ltd. nor Barry Smee accept any liability for any decisions or actions taken following the use of the reference material. Our liability is limited solely to the cost of the reference material.

Certified by

  
\_\_\_\_\_  
Duncan Sanderson, Certified Assayer of B.C.

Geochemist

  
\_\_\_\_\_  
Dr. Barry Smee, Ph.D., P. Geo.